



EmpowerMR[®]

Contrast Injection System



A unique system
delivering **optimal
flow** in the highest
Tesla fields

**BRACCO
INJENEERING**
Innovative Injection Technology



EmpowerMR®

A state-of-the-art contrast injection system

- Approved for 7T scanners
- No batteries to manage

The EmpowerMR® contrast injection system offers a full range of unique features to optimize safety while helping technologists achieve optimal images.





ADVANCING PATIENT SAFETY

- **Tilt lockout:** Helps to minimize risk of air embolism
- **Color-coded, illuminating handle:** Indicating that contrast or saline is flowing to the patient
- **GFR calculator:** Helps identify patients with compromised renal function

ENHANCING YOUR PRODUCTIVITY


- **Intuitive, user-friendly interface:** Easy operation
- **Compact, rotating head:** More operational flexibility
- **Dual initialize and replace syringe:** Increases productivity
- **Non-magnetic injector head construction:** No limits on proximity to scanner
- **Battery-free hydraulic drive system:** Continuous power supply avoids risk of battery failure



Fast input



Simple clear function

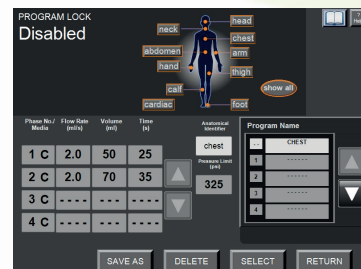

Remote screen



ORDERING INFORMATION



Code: 640471

EmpowerMR® CONTRAST INJECTION SYSTEM

Includes: Injector head with hydraulic controller, touch screen remote (Windows® operating system based) with desktop mount, electro-hydraulic interconnect, communication and power cables (includes pedestal mount for injector head).



Code: 017356

FAST-LOAD™ MR SYRINGE PACK, DCV

2 disposable 100 mL Fast-Load Syringes with 96" (y) low-pressure connecting tube with 2 x check valves. Contains 1 x large spike (for Saline bags), and 1 x small spike (for Contras Media bottles)

For more information on our complete contrast injector technologies portfolio, please visit our website <http://imaging.bracco.com>

CE 2797

© 2020 Bracco Injengineering SA. All rights reserved. P/N: TRCO-18872

Not all products are available in all global markets. Please refer to your local distributor or agent. Empower MR is a registered trademarks of Bracco Injengineering SA.

Bracco Injengineering SA reserves the right to modify the specifications and features described herein, or discontinue manufacture of the product described at any time without prior notice or obligation. Please contact your authorized Bracco Injengineering SA representative for the most updated information.

For additional information, please contact:

Bracco Injengineering SA
Avenue de Sévelin 46
1004 Lausanne
Switzerland

info.injengineering@bracco.com
<http://imaging.bracco.com>

BRACCO
INJENGINEERING
Innovative Injection Technology