

*Mirage*  
*Laguna*

**TPC**





# M4000 Hydraulic Patient Chair



## Articulating dual axis headrest

Headrest rotation allows for optimum positioning of the oral cavity without the need to reposition the chair



## Swing-away armrests

Help easy entry and exit from the chair





Seat rotates 30° to left or right from center to accommodate different working styles and to maximize access to intraoral x-ray system

## Mirage

## Luxury Design

TPC Mirage series dental chair have been sold in the USA for more than 10 years. It pays attention to the scientific, stable and high economy of product. Now we introduce this product to worldwide with its advantage simply, fashion and durable, you can enjoy the American Luxury dental chair with low cost.



### Linked Compensation

Seat inclination will be adjusted according to backrest movement, to keep patient in a best position.

## Safe and Reliable

1. Chair will rise automatically when encountering obstacles.
2. All functions will be temporarily turned off when using handpiece.
3. Any button on the foot control can cancel chair movement.

## 3 years warranty



Hydraulic motor system speed can be adjusted to ensure the patient remains comfortable as the chair changes position smoothly and silently.

Hydraulic motor provides smooth and quiet movement max. weight load 220kg.

# Airage

- Microfiber leather upholstery from USA, comfortable and durable.
- Seam or seamless are available.
- 8 Upholstery color selection.




Seam



Seamless



**DU-06**  
Desert Rose

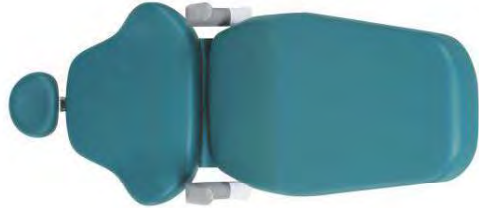

**DU-15**  
Grey Teal



**DU-18**  
Taupe



**DU-19**  
Greystone



**DU-100**  
Black



**DU-201**  
Atlantis



**DU-300**  
Cobalt




**DM02**  
Lake Green





## Mirage Doctor Stool / Assistant Stool

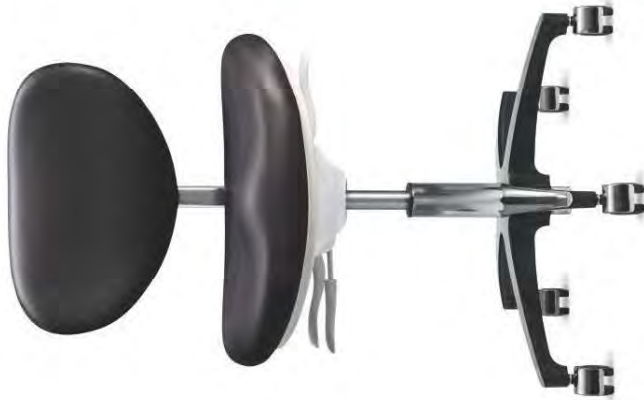

**DR3102**  
Doctor Stool

- Curved saddle seat and backrest can be adjusted to provide a comfortable sitting posture for long time operation.
- Cast aluminum base is solid and durable.
- 8 Upholstery color selection.



**AS3101**  
Assistant Stool (Optional)

- Stool arm can be rotated 360° easier for assistant to help doctor for operation.
- Foot pedal height is adjustable, fit for different sitting posture.
- Cast aluminum base with durable wheels ensures smooth movement.

*Mirage*

**MP4000(A)**  
Standard Type  
(Holder Type)



*Airage*

**MP4000(B)**  
Continental Type  
(Up Mounted Type)





# MP4000SW (Swing Mounted System)



MP4000SW(B)



MP4000SW(A)

Standard type instrument tray is available

Rotatable instrument holder arm pivots from under the seat, eliminating the pressure on patients caused by the arm of the instrument tray.



## TPC's luxurious and multifunctional dental units

- Durable tubing of air, water and handpiece, solenoid valve and 3 way syringe are from USA.
- Anti-suction system for handpiece tubing to prevent cross infection.



### Optional item:


- Equipped with tubing disinfection system.
- Patented design, tasteless, super sterilization function.



- 180° rotatable spittoon base.
- Flexible spittoon easy to clean.
- Electronic control of spittoon flushing time.
- The assistant position also comes with a full function control panel.

*Mirage*

## MP4000 Luxury Design

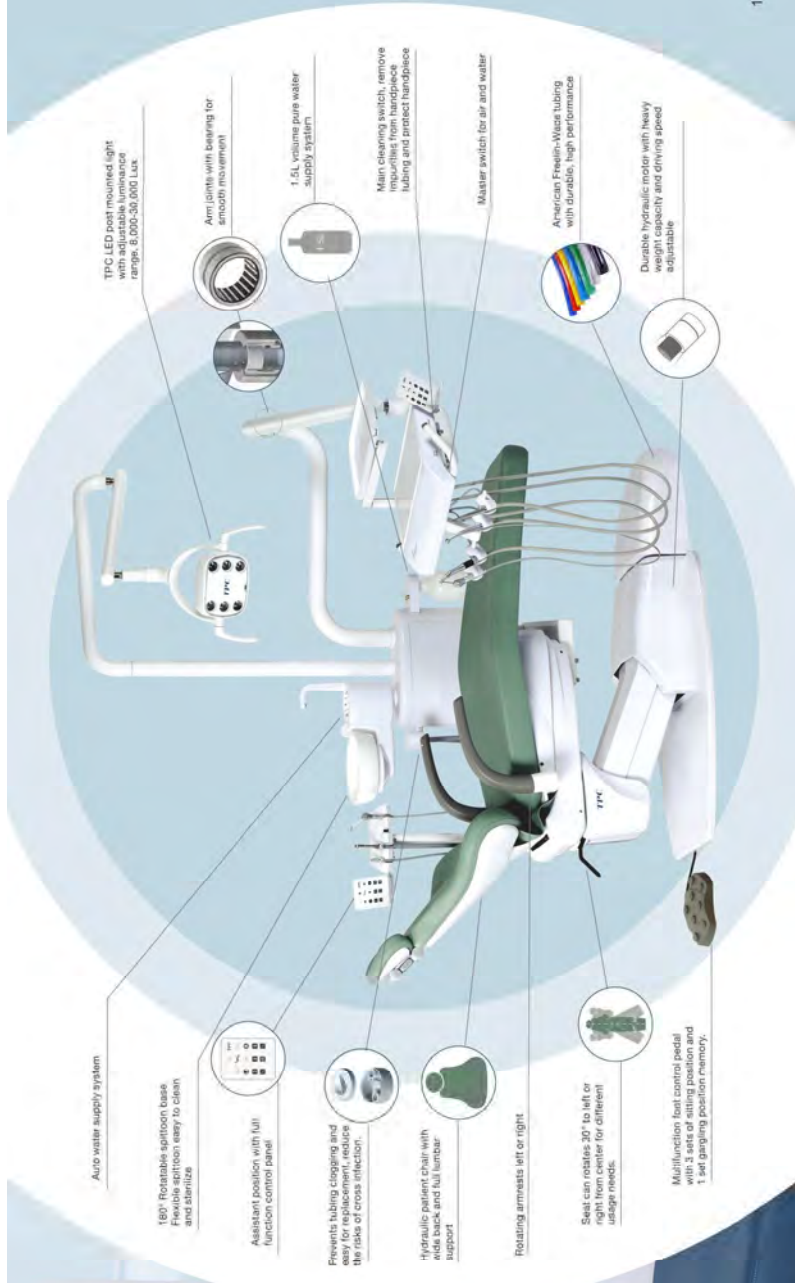

Auto water supply system

180° Rotatable spittoon base  
Flexible spittoon easy to clean  
and sterilize

Assistant position with Multi  
function control panel

Prevents tubing clogging and  
easy for replacement, reduces  
the risk of cross infection.

Hydraulic patient chair with  
back and full lumbar  
support

Rotating armrests left or right.

Seat can rotate 30° to left or  
right from center for different  
usage needs.

Multi-function foot control pedal  
with 3 sets of sitting position and  
1 set gangling position memory.

TPC LED post mounted light  
with adjustable luminance  
range: 8,000-30,000 Lux.

Arm joints with bearing for  
smooth movement

2.5L volume pure water  
supply system

Main cleaning switch, remove  
impurities from handpiece  
tubing and protect handpiece

Master switch for air and water

American Freeen-Wies tubing  
with durable, high performance

Durable hydraulic piston with heavy  
weight capacity and driving speed  
adjustable